



Revypeel

The lifting & bio-revitalising peeling

Trichloroacetic acid | Kojic acid | Mandelic acid

mbe

REVPYPEEL is a CLASS IIa medical device that performs skin cell renewal for face, neck and décolleté and body by providing deep skin stimulation without inducing any superficial irritation.

This chemical peel contains trichloroacetic acid, mandelic acid and kojic acid with an exfoliating action, specifically formulated to treat chrono- and photo-ageing, skin problems caused by acne and skin dyschromia, stretch-marks and scars.

Revypeel LOW

A light, superficial chemical peel that acts on the stratum corneum and stratum granulosum of the epidermis for treating mild acne lesions (comedones) and superficial scarring. This mix of acids provides gentle exfoliation with a depigmenting action to treat mild dyschromia. By removing the most damaged outer layers, cell renewal is promoted and wrinkles, surface scars and skin spots are reduced.

Composition:

Revypeel Low contains the following ingredients at different concentrations:

- Trichloroacetic acid (12%)
- Mandelic acid
- Kojic acid
- Excipients



Revypeel HIGH

A chemical peel with a high concentration of trichloroacetic acid, which acts on the entire thickness of the epidermis, specifically designed to treat hypertrophic cheloids, deep scars, stretch marks, moderate to severe acne lesions (papules, pustules, nodules and cysts) and related scarring. This mix of acids provides intense exfoliation with a depigmenting action to treat severe dyschromia. By removing the most damaged outer layers, cell renewal is promoted and wrinkles, medium-depth and surface scars, and skin spots are reduced.

Composition:

Revypeel High contains the following ingredients at different concentrations:

- Trichloroacetic acid (20%)
- Mandelic acid
- Kojic acid
- Excipients

Unlike classic biostimulation treatments, REVYPEEL does not need to be injected.

Active Ingredients

Mandelic Acid


Moisturising, Regenerating, Bactericidal, Nourishing Action

Trichloroacetic Acid (TCA)

Keratolytic, Lightening, Bactericidal, Stimulating Action

Kojic Acid

Brightening, Exfoliating, Moisturising Action

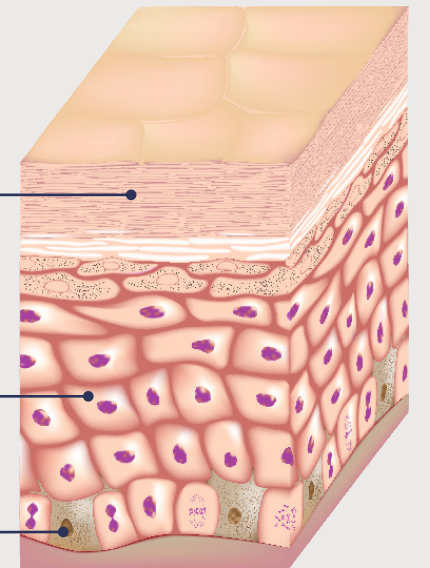

Features

Needle-free

No 'frost' effect

Suitable for all phototypes (I-VI)

Regenerative and revitalising stimulus for tissue self-repair and remodelling

Active in depth giving the skin a more compact, smooth and radiant appearance

Treatment to be performed only by medical professionals or under medical supervision

Indications

- **CHRONO-AGEING:**
It leaves the tissues firmer, smoother and brighter
- **PHOTO-AGEING:**
It helps remove spots and reduces skin discoloration
- **ACNE:**
It exerts a soothing anti-inflammatory action on acneic skins
- **DYSCHROMIAS:**
It lightens skin spots and minimizes dyschromia thanks to tyrosinase inhibition and the lightening of existing melanine
- **SCARS:**
It has a strong smoothing power on scar tissue
- **STRETCHMARKS:**
It evens out the and has a strong smoothing power

REVYSKIN Complementary Line

Pre-peeling lotion

Skin toning lotion. The special pool of acids contained in the lotion, evens the skin pH before proceeding with peeling treatment.

Main ingredients

ASCORBIC ACID (vit C.)

CITRIC ACID

LACTIC ACID

Neutralizing gel

It is an alkaline solution at pH9 designed to neutralize the action of REVYPEEL

Main ingredient

SODIUM CARBONATE

Post-peeling cream (SPF 50+)

Protective cream with wide spectrum uva and uvb filters, ideal for skin stressed from chemical peels. Prevents the formation of skin spots and reduces the damage of photoaging.

Main ingredients

PANTHENOL

HYALURONIC ACID

KARITÉ BUTTER

ARGAN OIL

VITAMIN E



PRO
TO
CO
L



PREPARATORY SOLUTION

Lotion formulated to standardise the skin PH. Its composition includes a cocktail of soft organic acids, that regulate the skin PH without attacking the epidermis.

1



APPLICATION OF REVYPEEL

After application it is necessary to deeply penetrate the active ingredients with a light manual massage

2



NEUTRALISING SOLUTION

To block the activity of the chemical peel it is necessary to buffer a specific solution.


3



SOOTHING CREAM

Scientifically studied cream with sunscreen, this cream is an excellent ally to protect, soothe and moisturize the epidermis stressed by chemical agents.

4

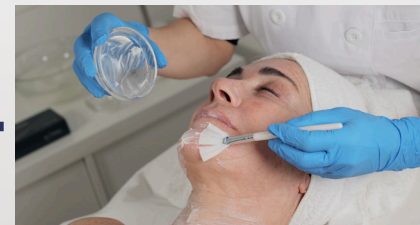

Used with a regular frequency, **REVYPEEL** is able to slow down the effects of intrinsic and extrinsic skin ageing, while obtaining results right after the first session.

- Contrasts chrono and photo-ageing
- Reduces fine wrinkles
- Revitalises dermis and epidermis
- Lightens and brightens the complexion
- Enhances hydration
- Anti acne both in the active phase and on scar outcomes
- Reduces stretch-marks
- Improves recent depressed scars
- Reduces the chloasma

Effects

- Bio-revitalizing
- Exfoliating
- Brightening
- Lightening
- Bactericidal

Synergies

KIND OF TREATMENT	WHEN TO USE REVYPEEL
INJECTABLES	BEFORE
PHOTO STIMULATION	BEFORE
NEEDLING	BEFORE
TOPICAL TREATMENTS	BEFORE
RADIOFREQUENCY / LASER	AFTER

Proven Results

Wrinkles & Photo-ageing

Revypeel LOW



Fine wrinkles

Revypeel LOW



Wrinkles

Revypeel LOW



Skin Spots

Revypeel LOW



Scars

Revypeel LOW



Chrono-ageing

Revypeel LOW



Acne

Revypeel LOW



Striae

Revypeel HIGH



mbc

mbmedicale.it

MBE is a G. TRADING COMPANY BRAND

Head Office: Via Castiglione 99, 40136 BOLOGNA

Warehouse: Via Bruno Tosarelli 356, 40055 Villanova di Castenaso (BO)

info@gtrading.it

Disclaimer: All contents of this material are for informational purposes only and provided by G TRADING COMPANY without warranties of any kind. Healthcare professionals are solely responsible for making their own independent evaluation as to the suitability of any product for any particular purpose and in accordance with country specific regulations. The availability of products and the indications mentioned in this material are subject to the regulatory requirements and product registration status in each country. Refer to your User Manual for your country specific indications. Products and technical specifications may change without notice. Please contact G TRADING COMPANY for more details.

© 2024 G TRADING COMPANY. This material contains registered and unregistered trademarks, trade-names, service marks and brand names of G TRADING COMPANY and its affiliates, including G TRADING COMPANY. All other trademarks are the property of their respective owners. All rights reserved.

